



**WAARDENBURG**  
MEDICAL

**Derungs**  
MEDICAL LIGHTING



TRIANGO 60

The compact one that has what it takes



# TRIANGO 60 – your new treatment & minor surgery light



**TRIANGO 60 C**  
ceiling version




**TRIANGO 60 W**  
wall mount



**TRIANGO 60 F**  
roller stand



**TRIANGO 60 F Battery**  
roller stand

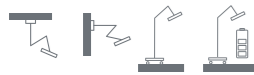

| TRIANGO 60                     | technical data                                    |
|--------------------------------|---|
| Fitted with                    | 1x LED-Modul                                      |
| Luminous power *               | 60 000 Lux / 1,0 m                                |
| Colour temperature *           | fixed value 4 700 K                               |
| TRIANGO 60-1 F                 | changeable value                                  |
| TRIANGO 60-3 F                 | 3 700 / 4 300 / 4 700 K                           |
| Lighting field 1,0 m*          | ø 18 cm   |
| Depth of illumination *        | 140 cm  |
| Colour rendering index *       | Ra > 94   |
| Operation equipment            | Electronic converter                              |
| Voltage                        | 100-240 V; 50/60 Hz                               |
| Luminaire body                 | Polyamide (PA)                                    |
| Power consumption              | < 30 W  |
| Usage                          | Touch panel on cardanic system                    |
| Technology                     | Single-pole secondary                             |
| Balance of articulation joints | Spring (weight compensation)                      |
| Class of protection            | I   |
| Protection rating              | IP 43 (illumination head)<br>IP 20 (UPSmedv Pack) |
| Standards met                  | EN 60601-1, EN 60601-2-41                         |

| Special features              |  |
|-------------------------------|--|
| <b>TRIANGO 60-1</b>           | <ul style="list-style-type: none"> <li>Continuously dimmable</li> </ul>                                    |
| <b>TRIANGO 60-3</b>           | <ul style="list-style-type: none"> <li>Continuously dimmable</li> <li>3 colour temperatures</li> </ul>     |
| <b>TRIANGO Fokus 60-1</b>     | <ul style="list-style-type: none"> <li>Continuously dimmable</li> <li>Focusable</li> </ul>                 |
| <b>TRIANGO 60-1 F Battery</b> | <ul style="list-style-type: none"> <li>Continuously dimmable</li> <li>UPSmed Pack</li> </ul>               |
| <b>Mounting</b>               | <ul style="list-style-type: none"> <li>Ceiling mount</li> <li>Wall mount</li> <li>Roller stand</li> </ul>  |
| <b>Cord length</b>            |  |
| <b>TRIANGO 60-1 C</b>         | Direct power supply  |
| <b>TRIANGO 60-1 W</b>         | Direct power supply  |
| <b>TRIANGO 60-1 F</b>         | 3 m  |
| <b>TRIANGO 60-1 F Battery</b> | 3 m  |
| <b>Weight (net)</b>           |  |
| <b>TRIANGO 60-1 C</b>         | 13 kg  |
| <b>TRIANGO 60-1 W</b>         | 12,5 kg  |
| <b>TRIANGO 60-1 F</b>         | 19 kg  |
| <b>TRIANGO 60-1 F Battery</b> | 31 kg  |
| <b>Accessories</b>            | <ul style="list-style-type: none"> <li>Sterilizable handle</li> <li>Disposable cover for handle</li> </ul> |

**Versatile**  
versions with focusable or colour changeable illumination field

**Homogeneous light field**  
thanks to unique diamond optic

**Easy cleaning**  
streamlined design and the removable handle

**Hygiene made simple**  
good design does without unnecessary details and attaches importance to closed building structures

**Can be used everywhere**  
TRIANGO 60 is suitable for the versatile daily routine in family practices and clinics as well as in gynecology, dermatology or veterinary medicine



**Ergonomic handling**  
easy alignment for right- and left-handed



**Intuitive operation**  
thanks to an intuitive touch panel on the cardanic, the illumination from 60 000 lux at 1 m can be infinitely adjusted

| Order number: TRIANGO 60-1    |   |
|-------------------------------|---|
| <b>TRIANGO 60-1 C</b>         | D16.250.000-00813741<br>D16.253.000-00814312 - USA  |
| <b>TRIANGO Fokus 60-1 C</b>   | D16.252.000-00814311<br>D16.255.000-00814314 - USA  |
| <b>TRIANGO 60-1 W</b>         | D16.256.000-00813742<br>D16.274.000-00814327 - USA  |
| <b>TRIANGO Fokus 60-1 W</b>   | D16.258.000-00814316<br>D16.276.000-00814329 - USA  |
| <b>TRIANGO 60-1 F</b>         | D16.262.000-00813744 - CH<br>D16.259.000-00813743 - Schuko<br>D16.271.000-00813747 - UK<br>D16.265.000-00813745 - USA |
| <b>TRIANGO Fokus 60-1 F</b>   | D16.264.000-00814320 - CH<br>D16.261.000-00814318 - Schuko<br>D16.273.000-00814326 - UK<br>D16.270.000-00814322 - USA |
| <b>TRIANGO 60-1 F Battery</b> | D16.279.000-00814482 - CH<br>D16.278.000-00814481 - Schuko<br>D16.282.000-00814485 - UK<br>D16.280.000-00814483 - USA |

| Order number: TRIANGO 60-3 |   |
|----------------------------|---|
| <b>TRIANGO 60-3 C</b>      | D16.251.000-00814310<br>D16.254.000-00814313 USA  |
| <b>TRIANGO 60-3 W</b>      | D16.257.000-00814315<br>D16.275.000-00814328 USA  |
| <b>TRIANGO 60-3 F</b>      | D16.263.000-00814319 - CH<br>D16.260.000-00814317 - Schuko<br>D16.272.000-00814325 - UK<br>D16.266.000-00814321 - USA |



\* -10% / +20% tolerance



## Headquarters

Derungs Licht AG  
Hofmattstrasse 12  
9200 Gossau  
Switzerland  
+41 71 388 11 66  
mailbox@derungs.swiss  
www.derungs.swiss




**WAARDENBURG**  
MEDICAL

Dorpsstraat 8a  
4181 BN Waardenburg

www.waardenburgmedical.com  
info@waardenburgmedical.com

+31(0) 418 651 078  
+31(0) 61066 2894

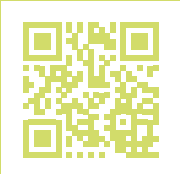

Derungs Licht AG is  
a company of the  
Waldmann Group

[www.derungs.swiss](http://www.derungs.swiss)



Register and benefit:  
+ 1 year warranty and  
exclusive offers

[www.derungs.swiss/product-registration](http://www.derungs.swiss/product-registration)